



The new frontier in diagnostic hysteroscopy

Code	MC.4SP-03004-001
Description	Kit composed of: - 1 Just-On - 1 Sterile non-woven gauze, 10x10cm - - 4 layers

Clinical reliability, ease of use and broad compatibility make JuSt-On an essential tool for hospitals, specialist centres and gynaecology practices.



Information:

 Zona Industriale
67062 Magliano dei Marsi (AQ)

 Ph.: +39 0863 517956

 info@assuteurope.com

 www.assuteurope.com

Assut Europe S.p.A.

Founded in Italy in 1991, Assut Europe S.p.A. is a leading manufacturer of surgical sutures and medical devices for a wide range of surgical specialties, including general, vascular, cardiovascular, thoracic and orthopaedic surgery.

Present in over 11 countries worldwide, it stands out for the quality of its products and a vision based on the excellence of "Made in Italy", with constant investment in innovation and research.




SCAN ME



©2026 Assut Europe S.p.A. All rights reserved

Just-on

Innovation that goes beyond blind biopsy

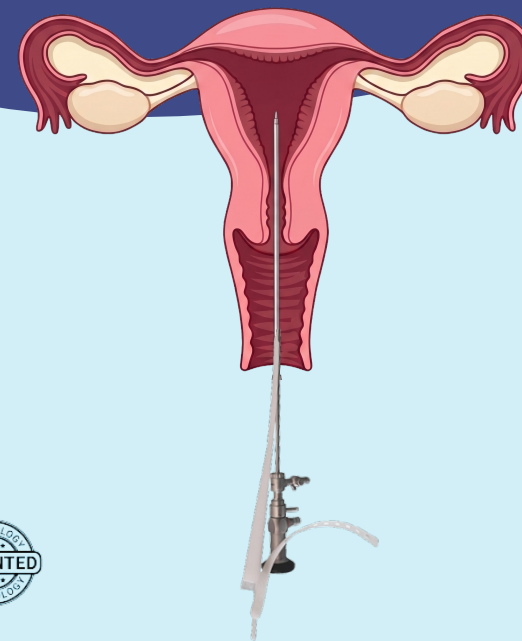

Turns diagnostic hysteroscopy into an operative procedure

Precision, safety and versatility to turn any diagnostic hysteroscope into an operative instrument

JuSt-On is an innovative and patented medical device designed to integrate seamlessly with standard diagnostic hysteroscopy systems. It allows professionals to perform histological sampling, targeted biopsies, and minimally invasive procedures without disrupting clinical workflow. Its design ensures maximum reliability and control at every stage of the procedure.

A more efficient, uninterrupted procedure



Ready to use **at any time**



Compatible with instruments up to **6 Ch**



Turns a diagnostic hysteroscope into an **operative** one



Can be inserted **without leaving the uterine cavity**



No interference with the distension flow



Drainage channel with an **integrated clamp**

A reliable device for daily clinical practice

1

Controlled stability

Anchoring system that ensures greater safety even during the most complex operating maneuvers.

3

Adaptability and smoothness

Compatible with 4 mm diagnostic hysteroscopes, reduced profile, and self-lubricating material for atraumatic access

2


Integrated working channel

Passage of instruments up to 6 Ch with constant flow for always optimal vision

4

Smart fluid management

Anti-reflux access port and drainage with integrated clamp for effective distension of the uterine cavity

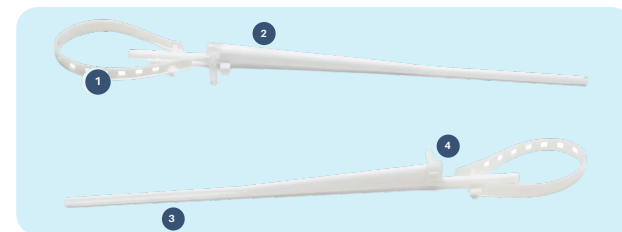

Step by Step

Fixation

With a light pressure, JuSt-On attaches to the hysteroscope, ensuring a stable and secure connection.

Positioning

The device is slid along the instrument until it reaches the desired operating position.

Instrument insertion

It is possible to insert the biopsy forceps or any other compatible instrument up to 6 Ch, according to clinical needs

Operative procedure

The sample is taken with full control maintained, without interrupting the hysteroscopic procedure.

